


HOMOGENIZE YOUR SAMPLES MORE PROTECTION

THE BIOBOX™

A Locking Containment Unit for the Bullet Blender®

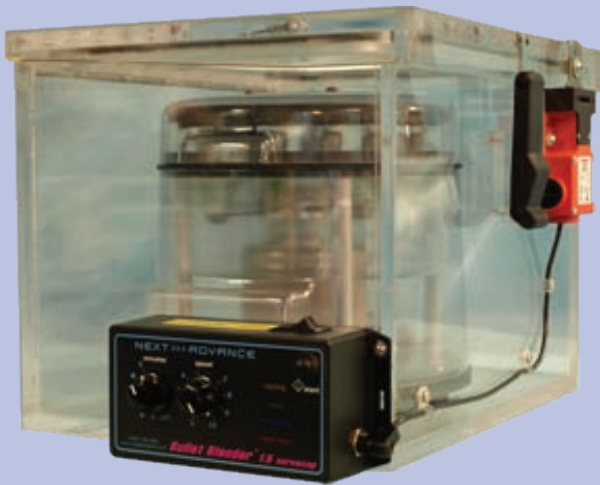

BIOBOX™ with BULLET BLENDER® inside

Now you can get an additional layer of protection when you homogenize pathogens or infected samples. The BIOBOX™ is a protective containment system that houses the BULLET BLENDER® homogenizer. The BULLET BLENDER® can homogenize tissue and cell culture for the isolation of a variety of analytes. The BIOBOX™/BULLET BLENDER® combination works with standard screw cap microcentrifuge tubes and retains the ease of use and high value that the BULLET BLENDER® is known for.

Features

- Clear housing so you can detect any leaks or spills caused by tube failure
- Microprocessor-controlled electronics integrate BIOBOX™ and BULLET BLENDER® into one easy to control unit
- Interlocks ensure that the BIOBOX™ is properly sealed
- Negative Pressure so that contaminants remain inside of the containment unit in the unlikely event of tube failure
- Control panel is safely and conveniently placed on the outside of the containment unit
- Fits inside most biosafety cabinets for additional protection

NEXT >>> ADVANCE

INNOVATIVE PRODUCTS • BETTER WORKFLOW
ACCURATE RESULTS • www.nextadvance.com

24 Prospect Avenue
Averill Park, NY 12018 USA
Telephone: 800-738-1681
Fax: 518-674-0189
info@nextadvance.com